



Domus 3

Alternating Pressure Redistribution System

- Suitable for patients with medium risk of pressure ulcers
- Ventilated low air loss for patients' utmost comfort
- Effective alternating system with static therapy mode in help with patient repositioning
- Low pressure alarm feature increase patient safety





Domus 3

Alternating Pressure Redistribution System



**Ventilated
Low Air Loss**
Wicks moisture
away to keep the
patient comfortable



Foam Base Pocket
For inserting additional
foam available on 5"
mattress
(Foam is not included)



**Zipper-all-Around
Cover**
Reduce ingress of
unwanted fluid and
simply one step to
remove for cleaning



CPR Valve
Quickly deflate the mattress
during emergencies



Nylon/PU Cover
Vapor permeable
and water resistant
with basic flame
retardant level



Quick Connector
Once with cap on,
instantly switch to
static air mattress for
transportation need



Low Pressure Alarm
Feature increase
patient safety



Biocompatibility



Vapour Permeable




Water Resistance



Flame Resistance

Pressure Relief During Alternating

Pressure mapping reveals that in alternating mode, bladder A & B take turns to inflate and deflate can significantly relieve the body pressure

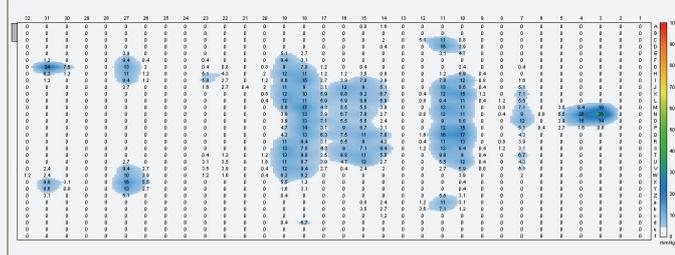

2 Modes to Fulfill Individual Needs

- Alternating
- Static



Specification

Model	Domus 3		
Pump	Dimension	11.4" x 7.3" x 5" / 29 x 18.5 x 12.6 cm	
	Weight	5.1 lb / 2.2 Kg	
	Case Material	Flame retardant ABS	
	Supply Voltage	AC 100-120V 60Hz, 0.16A	
	Operating Cycle	8 min	
Mattress	Dimension	78.7" x 35.4" x 5" w/ H3" pocket	78.7" x 35.4" x 8"
		200 x 90 x 12.7 cm w/ H7.6 cm Pocket	200 x 90 x 20.3 cm
	Cell Height	5" x 17 cells	8" x 20 cells
	Weight	12.5 lb / 5.7 Kg	15.2 lb / 6.9 Kg
	Material	Cover: Nylon/PU; Cell: Nylon/PVC; Base: Nylon/PVC	
Maximum Patient Weight	440 lb / 200 Kg		
Flame Retardant Standards	CAL 117		
HCPCS code	E0372	E0277	

Specifications are subject to change for improvement without notice

Ordering Information (For North America Only)

Model	5" Overlay + 3" Foam pocket	8" Replacement	
System	With Quilted Cover	PA01401	PA01501
	With Non-Quilted Cover	PA01402	PA01502
Mattress	With Quilted Cover	PM08401	PM09405
	With Non-Quilted Cover	PM08402	PM09406
Pump	PA01001		