



# Domus 4

## Alternating Pressure Redistribution System with Low Air Loss

- Suitable for patients with medium to high risk of pressure ulcers
- Ventilated low air loss for patients' utmost comfort
- Cell-on-cell bladder to avoid bottoming out
- Visual and audible alarm enhances patients safety and security





# Domus 4

## Alternating Pressure Redistribution System with Low Air Loss



**Ventilated Low Air Loss**  
Wicks moisture away to keep the patient comfortable



**Cell-on-Cell Bladder**  
The bottom cells acts as basic pressure support to avoid bottoming out in case of power failure



**Zipper-all-Around Cover**  
Reduce ingress of unwanted fluid and simply one step to remove for cleaning



**CPR Valve**  
Quickly deflate the mattress during emergencies



**Nylon/PU Cover**  
Vapor permeable and water resistant with basic flame retardant level



Biocompatibility



Vapour Permeable



Water Resistance



Flame Resistance



**Quick Connector**  
Once with cap on, instantly switch to static air mattress for transportation need




**Visual and Audible Alarm**

- Low Pressure
- Power Failure
- Tech Support


## Pressure Relief During Alternating

Pressure mapping reveals that in alternating mode, bladder A & B take turns to inflate and deflate can significantly relieve the body pressure

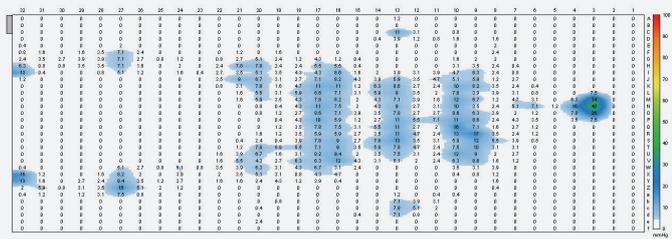

### 5 Modes to Fulfill Individual Needs

- Max Firm
- Alternating
- Continuous Low Pressure (CLP)
- Seat Inflation w/ Alternating
- Seat Inflation w/ CLP

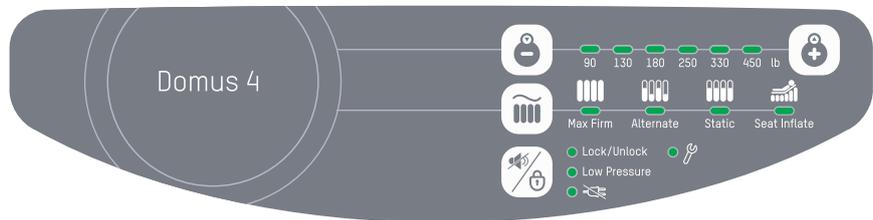

## Specification

<b>Model</b>	<b>Domus 4</b>	
<b>Pump</b>	Dimension	11.4" x 7.3" x 5.0" / 29 x 18.5 x 12.6 cm
	Weight	5.1 lb / 2.3 Kg
	Case Material	Flame retardant ABS
	Supply Voltage	AC 100-120V 60 Hz, 0.17A
	Operating Cycle	10 mins
<b>Mattress</b>	Dimension	8" Replacement
	Dimension	78.7" x 35.4" x 8" / 200 x 90 x 20.3 cm
	Cell Height	8" x 20 cells
	Weight	15.2 lb / 6.9 Kg
	Material	Cover: Nylon/PU; Cell: Nylon/PU; Base: Nylon/PVC
	Maximum Patient Weight	450 lb / 200 Kg
	Flame Retardant Standards	CAL 117
<b>HCPCS code</b>	<b>E0277</b>	

*Specifications are subject to change for improvement without notice*

## Ordering Information (For North America Only)

<b>Model</b>	<b>Domus 4</b>	
<b>System</b>	With Quilted Cover	PD02501
	With Non-Quilted Cover	PD02502
<b>Mattress</b>	With Quilted Cover	PM09403
	With Non-Quilted Cover	PM09404
<b>Pump</b>	PD02001	

## HEADQUARTERS

### Wellell Inc.

No. 9, Minsheng St., Tucheng Dist.,  
New Taipei City 236044, Taiwan  
T +886 2 2268 5568  
F +886 2 2268 6525  
marketing@wellell.com

## PAN-ASIA

### CHINA

Wellell (Kunshan)Co., Ltd.  
No. 1368, ZiZhu Road, Kunshan,  
Jiang Su Province, China  
T +86 512 8617 1660 Ext 1150  
F +86 512 8617 8663  
customer.cn@wellell.com

### THAILAND

Apex Medical (Thailand)  
No.111/152 Soi Pho Kao Yak 19,  
Nawamin, Beung Kum,  
Bangkok 10240,  
Thailand.  
T +66 02 018 2140  
sales.thai@wellell.com

## PAN-EUROPE

### SPAIN

Apex Medical S.L.  
Elcano 9, 6ª planta 48008  
Bilbao. Vizcaya. Spain  
T +34 94 4706408  
F +34 94 4706409  
info@apexmedical.es

### UK

Wellell UK Limited  
Unit 33, Great Western Business Park  
McKenzie Way, Worcester ,WR4 9GN  
United Kingdom  
T +44 (0) 190 5774 695  
sales.uk@wellell.com

### FRANCE

Wellell France S.A.S  
4 Boulevard de la Chanterie  
49124 Saint-Barthelemy d' Anjou  
France  
T +33 2 72 79 27 13  
F +33 2 85 73 07 67  
info.france@wellell.com

### GERMANY

SLK Vertriebsgesellschaft mbH  
Am Herdicksbach 18, 45731 Waltrop  
Germany  
T +49 (0) 231-925360-0  
F +49 (0) 231-925360-29  
info@slk-gmbh.de  
www.slk-gmbh.de

## AMERICAS

### USA

Wellell America Corp.  
927 Mariner Street,  
Brea, CA 92821, USA  
T +1 714 671 3818  
F +1 714 494 8173  
sales.usa@wellell.com

Distributed by